

What you need to do:

The detection efficiency calibration of real germanium detector arrays!



- Install GRSISort
- Fit gamma-ray photopeaks with GRSISort
- Edit and run the RootEffi script
- Submit your fitted peaks, your edited RootEffi script and your efficiency curve ☺

What you need to know:

- How semiconductor detectors work
- How to determine the detection efficiency
- How to find nuclear information
- How to compile a code
- How to fit peaks with GRSISort
- How to run a script in ROOT and GRSISort

*GRIFFIN,
TRIUMF,
Vancouver,
Canada*

64 crystals

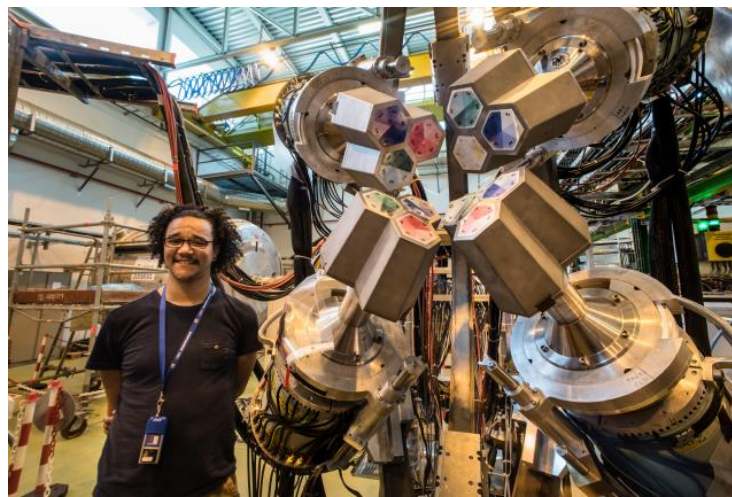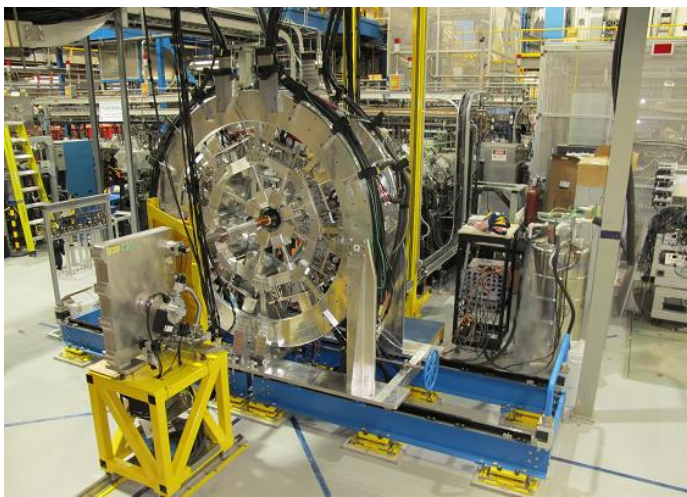


*MINIBALL,
CERN,
Geneva,
Switzerland*

24 crystals

*Soccer Ball,
iThemba LABS,
Cape Town,
South Africa*

52 crystals



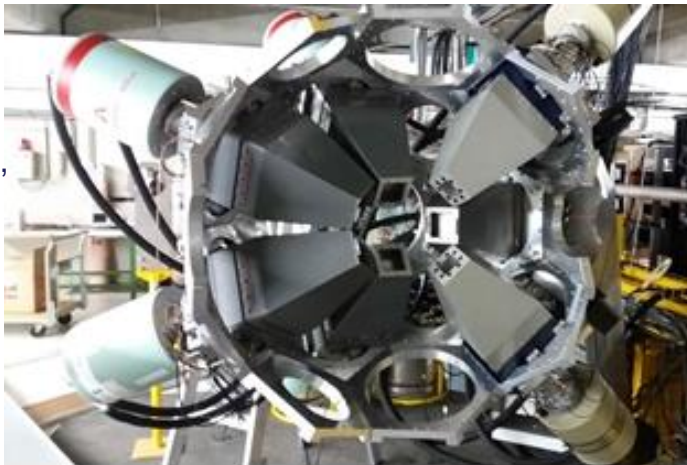


*TIGRESS,
TRIUMF,
Vancouver,
Canada*

56 crystals

GRIFFIN,
TRIUMF,
Vancouver,
Canada

64 crystals
(16 clovers,
not segmented)

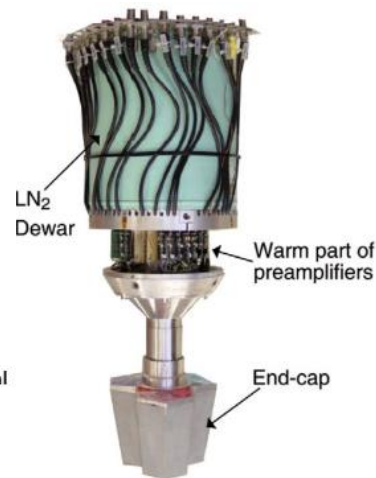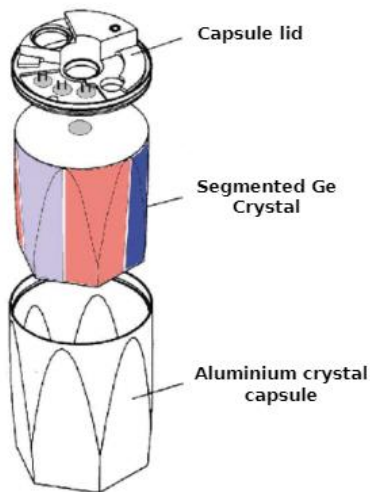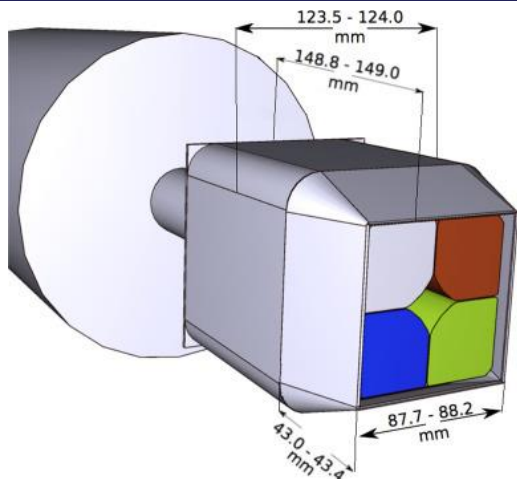



MINIBALL,
CERN,
Geneva,
Switzerland

24 crystals
(8 triple clusters,
6-fold
segmented)

Closed-end coaxial right cylinders ☺

Soccer Ball,
iThemba LABS,
Cape Town,
South Africa




52 crystals
(13 clovers,
4 segmented and
9 not segmented)

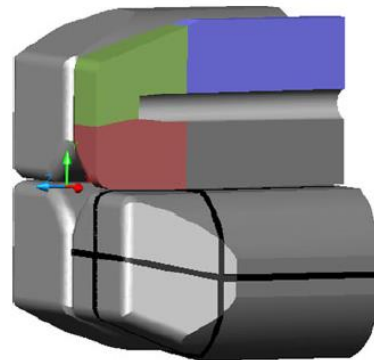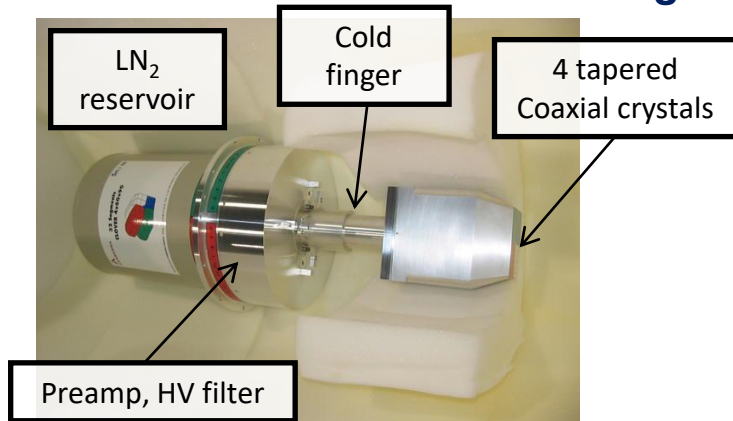

TIGRESS,
TRIUMF,
Vancouver,
Canada

56 crystals
(14 clovers,
8-fold
segmented)



- An external field creates an area depleted of free charge carriers.
- Radiation interacts with the crystal and produces electron-hole pairs.
- The electrons and holes drift towards electrodes and the electric pulse is amplified and processed.

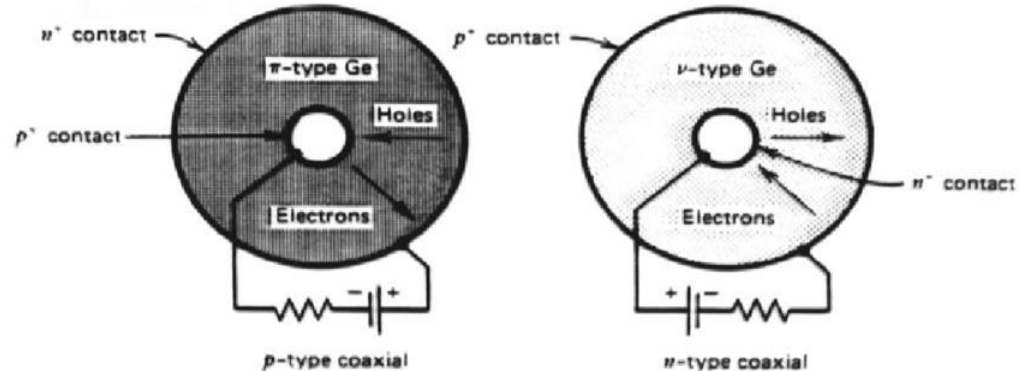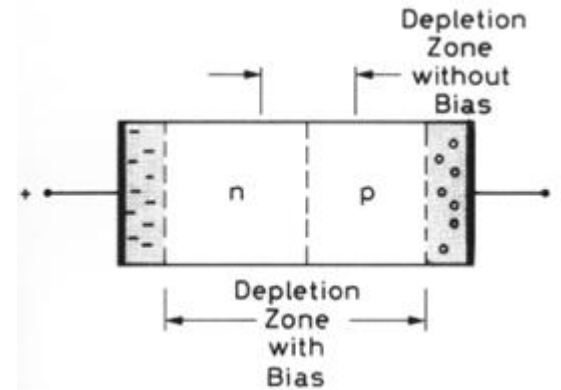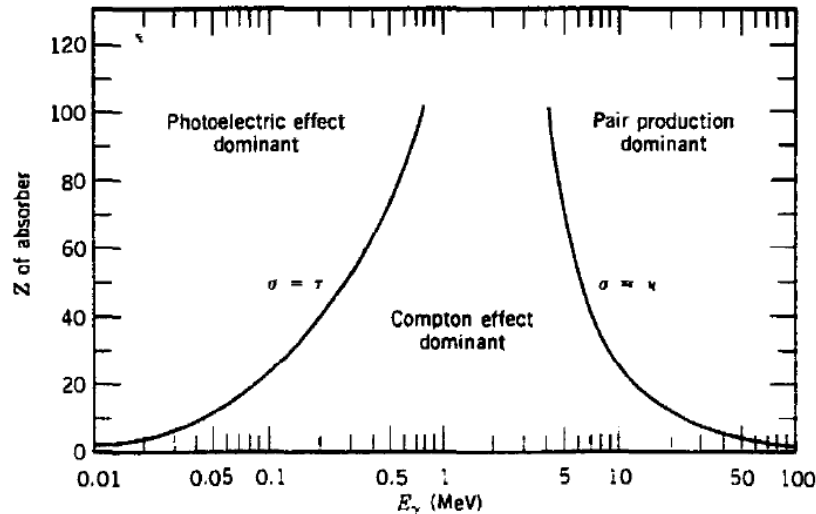


K.S. Krane, *Introductory Nuclear Physics* (1988).

G.F. Knoll, *Radiation Detection and Measurement* (1989).

W.R. Leo, *Techniques for Nuclear and Particle Physics Experiments* (1994).

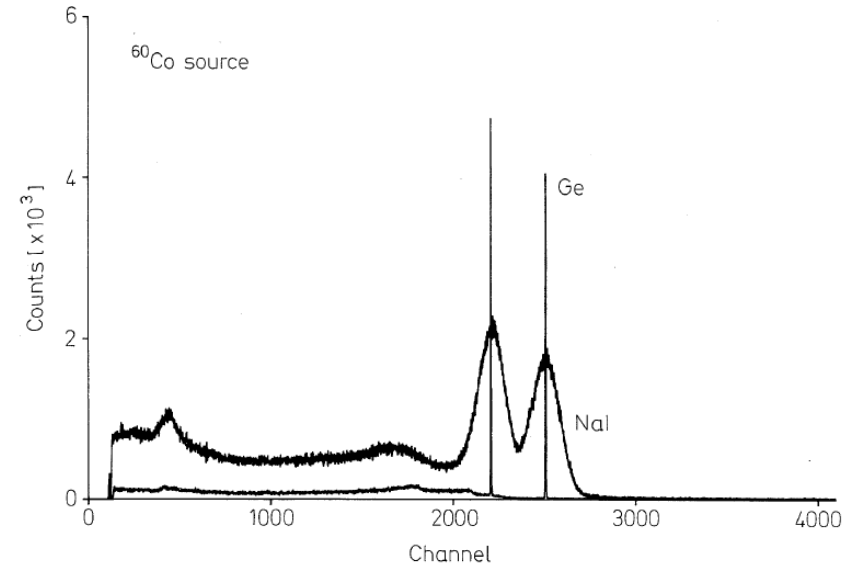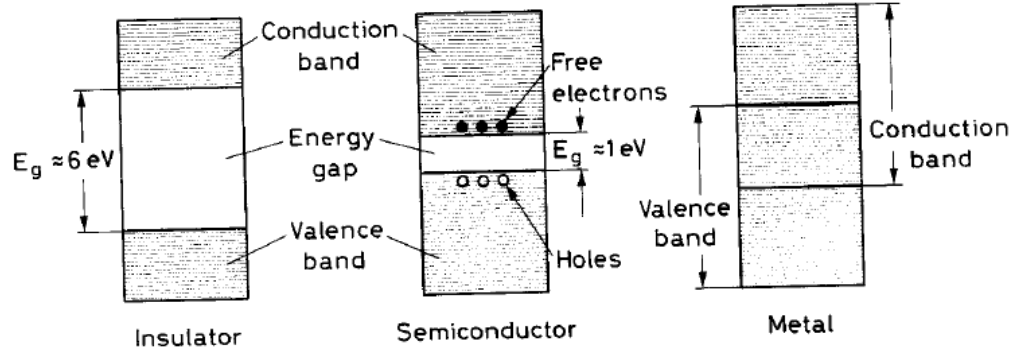

Material Property	NaI(Tl)	Germanium	Silicon
Type	Scintillator	Semiconductor	Semiconductor
Cooling	No	$\sim 95 \text{ K}$	No
Density (g/cm^3)	3.7	5.3	2.3
Band gap energy (eV)	20	2.9	3.8
Energy resolution at 1332 keV (keV)	13	1.9	

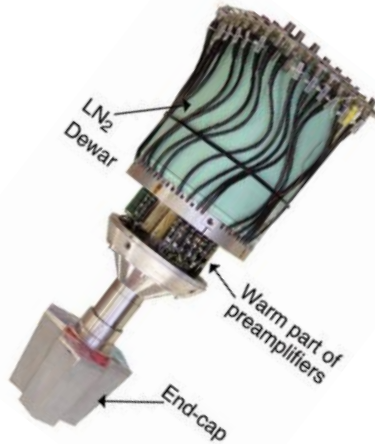

$$\text{Absolute Efficiency } \epsilon_{\gamma} = \frac{\text{Number of } \gamma\text{-rays detected}}{\text{Number of } \gamma\text{-rays emitted}}$$

$$= \frac{N_{\gamma, \text{detected}}}{I_{\gamma} \cdot A \cdot t}$$

← Area of the photopeak
 ← Data collection time
 ← Activity of the source
 ← Absolute intensity of the gamma ray

$$\epsilon_{\gamma}(E) = 10^{p_0 + p_1 \log(E) + p_2 \log^2(E) + p_3 / E^2}$$


 $z^A X$

Gamma(γ)decay

$$\text{Relative Efficiency } \epsilon_{\gamma} = \frac{\frac{N_{\gamma 1, \text{detected}}}{I_{\gamma 1} \cdot A \cdot t}}{\frac{N_{\gamma 2, \text{detected}}}{I_{\gamma 2} \cdot A \cdot t}} = \frac{N_{\gamma 1, \text{detected}} / I_{\gamma 1}}{N_{\gamma 2, \text{detected}} / I_{\gamma 2}}$$

- Lund/LBNL Nuclear Data Search: <http://nucleardata.nuclear.lu.se/toi/>
- National Nuclear Data Center: <https://www.nndc.bnl.gov/>

The Lund/LBNL Nuclear Data Search

Version 2.0, February 1999

S.Y.F. Chu¹, L.P. Ekström^{1,2} and R.B. Firestone¹

¹ LBNL, Berkeley, USA

² Department of Physics, Lund University, Sweden



WWW Table of Radioactive Isotopes

[Radiation search](#)

[Nuclide search](#)

[Atomic data](#) (X-rays and Auger electrons, very preliminary!)

[Periodic chart interface to the nuclides](#)

[Summary drawings for A=1-277](#) (PDF)

[Nuclear charts](#) (PDF, 333 kbyte)


[Database status](#)



Table of Isotopes (ToI)

[About this service](#)

[ToI home page](#)

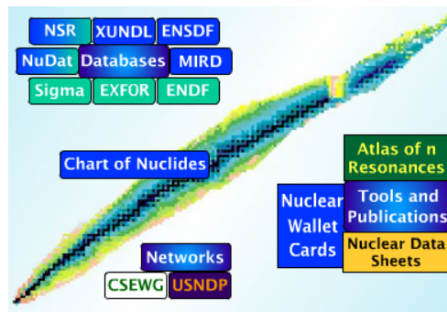

Tweets by @NNDC_BNL

National Nuclear Data Center

@NNDC_BNL

Meet the NNDC:

Boris Pritychenko is a experienced nuclear experimentalist and data compiler. He currently manages the EXFOR and NSR databases, and is also the Editor-in-Chief of Atomic Data and Nuclear Data Tables.

Embed

View on Twitter

Main Structure & Decay Reactions Bibliography Networks & Links Publications

AMDC Atomic Mass Data Center, *Q-value Calculator*

Covariances of Neutron Reactions

ENSDF Evaluated Nuclear Structure Data File

NSR Nuclear Science References

NuDat Nuclear structure & decay Data

Atlas of Neutron Resonances Parameters & thermal values

CSEWG Cross Section Evaluation Working Group

IRDF IRDF International Reactor Dosimetry and Fusion File

Nuclear Data Sheets Nuclear structure & decay data journal, *Special Issues* on reaction data

USNDP U.S. Nuclear Data Program

CapGam Thermal Neutron Capture γ -rays

EXFOR Nuclear reaction experimental data

MIRD Medical Internal Radiation Dose

Nuclear Wallet Cards Ground & isomeric states properties,

USNDP/CSEWG GForge Collaboration Server

Chart of Nuclides Basic properties of atomic nuclei

ENDF Evaluated Nuclear (reaction) Data File, *Sigma*

NDWG Nuclear Data Working Group

NucRates MACS & Astrophysical reaction rates

XUNDL Experimental Un-evaluated Nuclear Data List

- Installation instructions at <https://github.com/UWCNuclear/UbuntuSetUp>
- Open the file with “ grsisort -l Eu152_**ARRAYNAME**.root ”
- See the list of histograms in the file with “.ls”
- Draw histogram with “ gammaSingles->Draw() ”
- Display rough peak energies by pressing “s”
- Click and drag on x-axis to zoom in, zoom out by pressing “o”
- Click and drag on y-axis to zoom in, right click and click “UnZoom” to zoom out
- Click and click on spectrum to set the fitting boundaries around your peak
- Fit by pressing “f”
- Remove all markers by pressing “n”
- “.q” to quit
- Save screenshots of your fitted peaks and the number of counts (sum) in your peaks.
- More tools at <https://github.com/GRIFFINCollaboration/GRSISort/wiki/Interactive-Analysis>

- Download the data files with “ git clone <https://github.com/UWCNuclear/DetectorEfficiency.git> ”
- Download with “ git clone <https://github.com/UWCNuclear/RootEffi.git> ”
- In gedit, edit RootEffi.C with the peak areas and uncertainties obtained from your data file.
- Whenever you edit a script, you should save it and close ROOT/GRSISort to run it again.
- To run RootEffi, type in the command line “ grsisort -l RootEffi.C ”
or “ grsisort -l ” and then “ .x RootEffi.C ”
- Save your edited script and your new efficiency curve
- In one file, submit your fitted peaks, your edited RootEffi script, and your efficiency curve ☺